

Synergy ElectroSurgery System

High-Frequency Bipolar Ablation Console Quick Guide (AR-S9800)



Setup



1

Place the console on a flat, dry surface.



2

Connect the power cord to the console.



3

Connect the power cord to the outlet.



4

Power on the console.




5

Attach the probe via connector to the console.



6

Attach the footswitch on rear of the console.

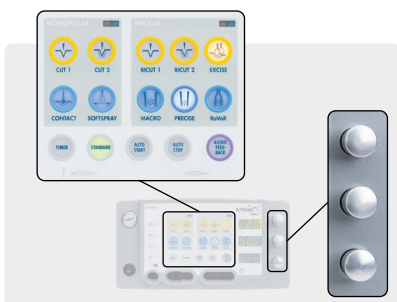

7

Change settings using the buttons and knobs on the front of the console.



Troubleshooting: Faults and Solutions

Fault	Solution
Unit does not power on	1. Check the main power plug and wall receptacle. 2. Call Arthrex Tech Support at 1-888-420-9393.
Accessory FlexTip RF probe does not work	1. Disconnect and reconnect the probe. 2. Turn the console off. Wait 30 seconds, then power it back on. 3. If applicable, disconnect and reconnect the footswitch. If the footswitch is not functioning, replace it. 4. Open new FlexTip RF probe.
'Err' message	1. Disconnect and reconnect the probe. 2. Restart the console. 3. Refer to the Synergy ElectroSurgery Console DFU (Ref 899800-E0) for error meaning, cause, and resolution.

To return or repair the device, contact Arthrex Product Surveillance at 1-866-267-9138.



Footswitch Activation	
1	Cut
2	Coag (Coagulation)