


PSI Osteotomy Case Upload Guide

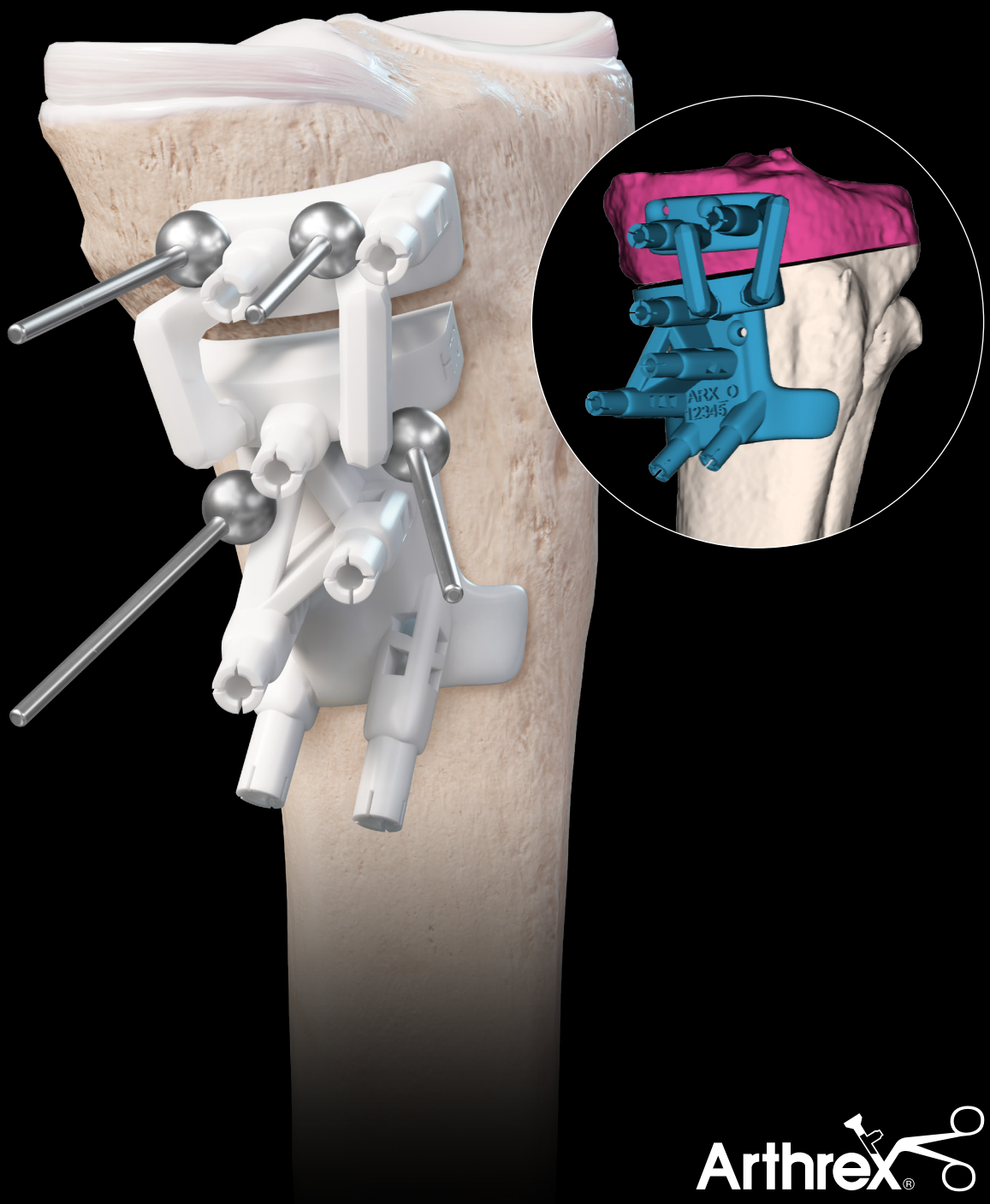


3D-Side Introduction

Arthrex partners with 3D-Side for case management and patient-specific instrumentation (PSI) for osteotomies. All healthcare providers and others requesting PSI instrumentation for upcoming osteotomy cases will use a secure platform hosted by 3D-Side for CT uploads and development of PSI instrumentation.

04	Website Account Creation or Login
04	New Users
05	Two-Factor Authentication
05	Email Authentication
06	Google Authentication App
06	Website Timeline View
07	Case Creation

Website Account Creation or Login

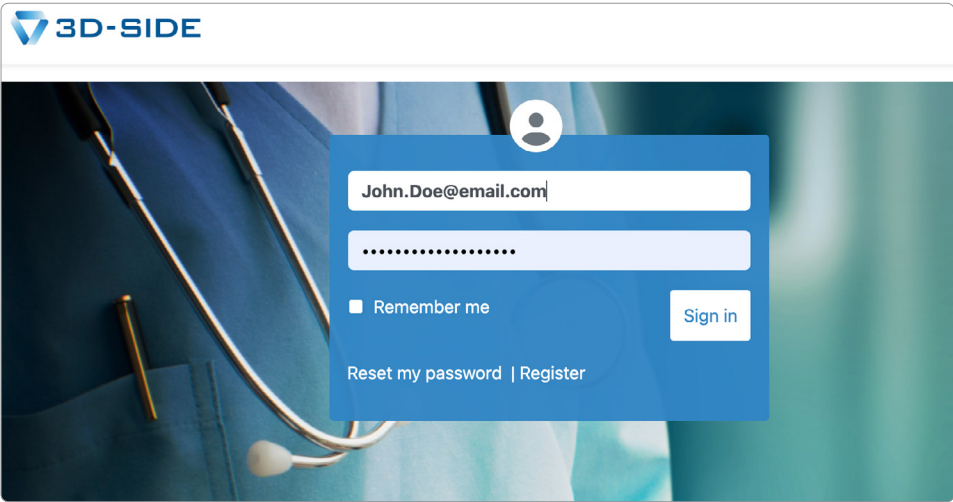

1

Navigate to <https://3dside.3d-customize.com/login>

| **Note:** Returning users will enter login information and skip to the “Timeline View” step (see page 6).

New Users

1

Click “Register” or navigate to <https://3dside.3d-customize.com/register> to set up an account.


Fill out basic information **(a)**:

- › Enter and confirm email
- › Create password
- › Select role, either “Surgeon” or “Sales Representative” **(b)**

2

Accept the data processing policy and click “Next” to fill out address information **(a)**.

Two-Factor Authentication

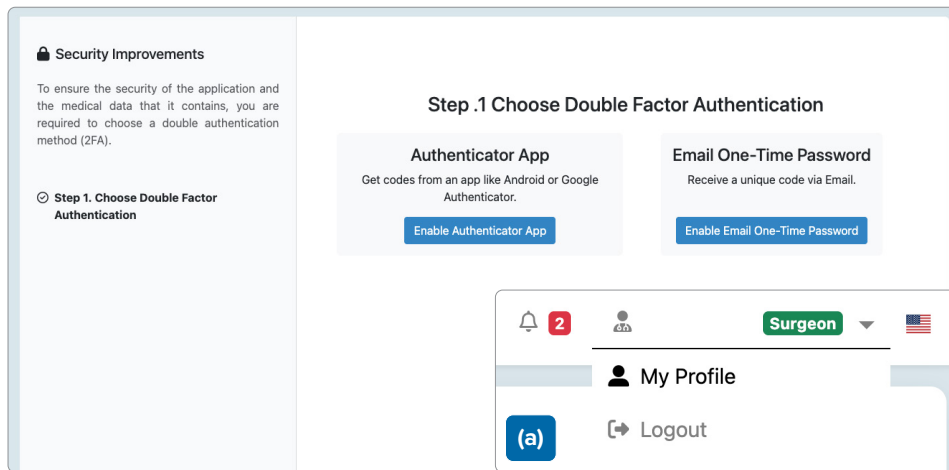


To comply with GDPR and HIPAA recommendations and reinforce the security of this platform, two-factor authentication is enabled.

- › New users will choose their preferred authentication method, either a six-digit code via email or the Google Authentication App.


Users can change their authentication method in their profile management page under “My Profile” at the top right of the page, under the user’s name **(a)**.

Email Authentication

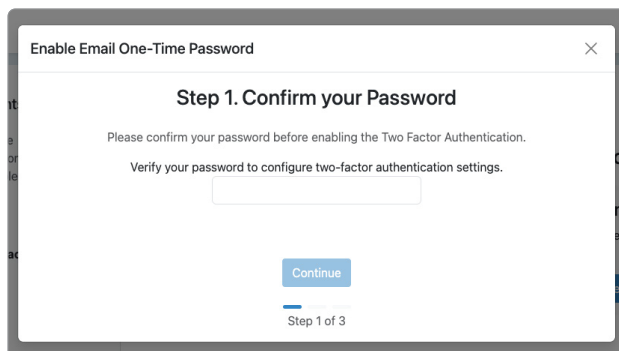


1

For email authentication, click “Enable Email One-Time Password,” then enter the password sent to the provided email address.

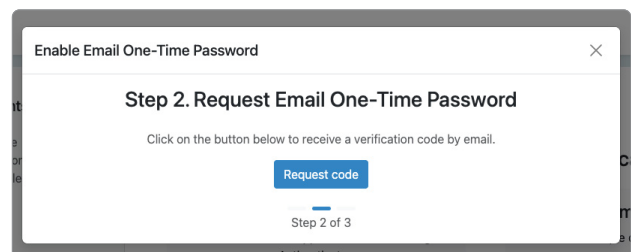

2


Click on “Request code.”

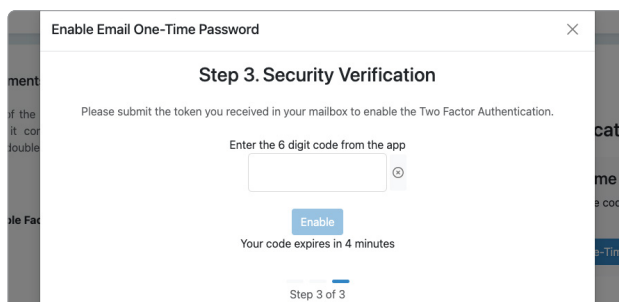

3


Enter the code sent to the provided email address.

Google Authentication App

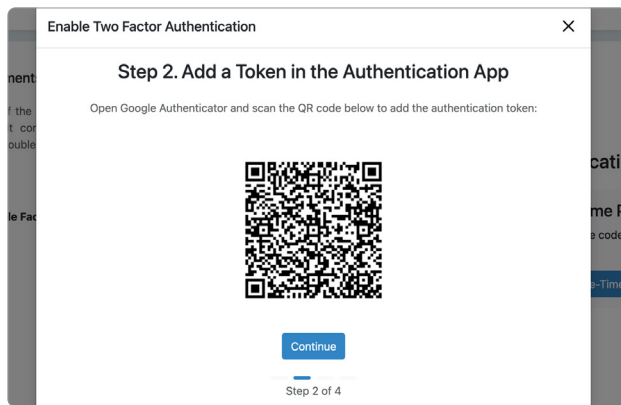


1

To use Google two-factor authentication, add a token in the authentication app by scanning the QR code displayed on the webpage then click “Continue.”

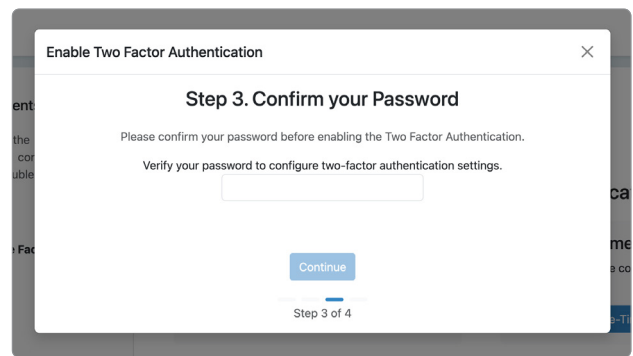

2

Confirm your password.

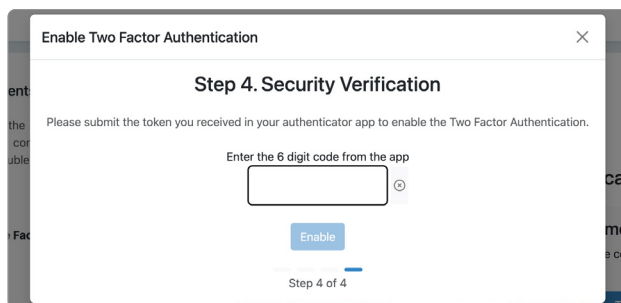


3


Enter the six-digit code found in the app into the space provided and click “Enable.”

Note: The code refreshes every 30 seconds.

Optional: If you have entered an incorrect code or need to switch accounts, click “Cancel Authentication Process” and restart the process.

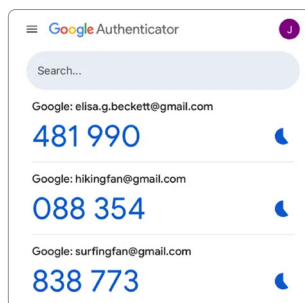


Website Timeline View

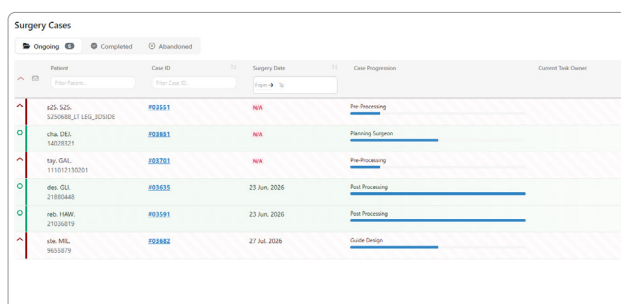

1

Bookmark <https://3dside.3d-customize.com/> to easily navigate to Timeline View, which is the standard landing page for users after creating an account.

- › Surgeon users will see their upcoming cases and due dates, with high-priority tasks clearly marked.
- › Arthrex sales representatives will see all currently pending cases for all surgeons.

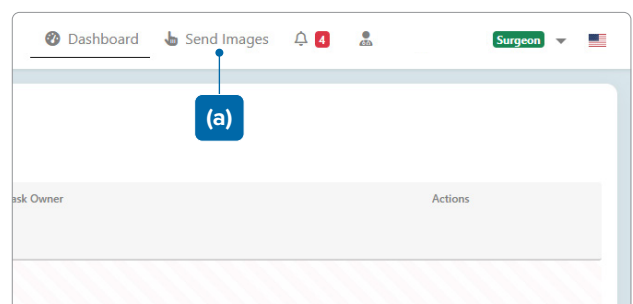

2

To create a new case and upload a new CT scan, click “Send Images” at the top right of the screen **(a)**.

Case Creation

3D-SIDE

Dashboard Send Images

Surgeon Information **Completed**

Surgeon First Name * Kevin ✓

Surgeon Last Name * Davis ✓

Surgeon Email * Kevin.Davis@Arthrex.com ✓

Surgery Information **Not completed**

Do You Know The Date Of Surgery?
☐ Yes ☒ No

Hospital *
 Type the first 3 letters of your Hospital

Pathology *
 Select a pathology

Additional Information *
 Please provide a description on the surgery and the specifications for the case

Please upload your medical images and complete the form.

Upload Your Images
 Complete the Form
 Create Series

Finish

Scan Protocol for Bone Tumor
 Scan Protocol for Cranioplasty
 Scan Protocol for Osteotomy

Drag and drop your images here. You can drop any files, zip or folders.

Add files

1

Navigate to <https://3dside.3d-customize.com/upload>

› For proper CT scan protocol, visit <https://rep.arthrex.com/resources/mis/DOC1-002268-en-US?t=1786546196649>

Surgeon Information **Completed**

Surgeon First Name * John ✓

Surgeon Last Name * Doe ✓

Surgeon Email * John.Doe@email.com ✓

Surgery Information **Not completed**

Do You Know The Date Of Surgery?
☐ Yes ☒ No

Hospital *
 Type the first 3 letters of your Hospital

2

If using a preregistered account, surgeon information will be prefilled. Otherwise, fill out surgeon information including the surgery details.

Note: The surgery date must be at least 4 weeks in the future to allow time for the PSI guide to be created and shipped.

Pathology * (a)
 Osteotomy

Surgical Side * (b)
 Right

Location * (c)
 Select location
 PAO
 PFO
 DFO
 HTO
 DTO
 Other
 be considered for the case

3

Under “Pathology,” the user must select “Osteotomy” (a), then “Surgical Side” (left or right) for the scan to be uploaded (b).

Select “Location” to choose a procedure type (ie, DFO or HTO) (c).

Surgical Side *
 Right

Location *
 HTO

Additional Information *

- Angles Of Correction
- Level Of Osteotomy
- HW Reconstruction Requested
- Surgical Approach

Please provide medical recommendation to be considered for the case


Fixation
 Reference your fixation if you have any preference (a)

4

Add case details such as desired correction angle, concomitant procedures, and/or patient history in the “Additional Information” box.


If a specific fixation method is desired (ie, the ContourLock™ plate, iBalance® system, or 4-hole technique), add the fixation method in the “Fixation” box (a).

Note: If left blank, the 3D-Side design team will suggest an appropriate fixation method.

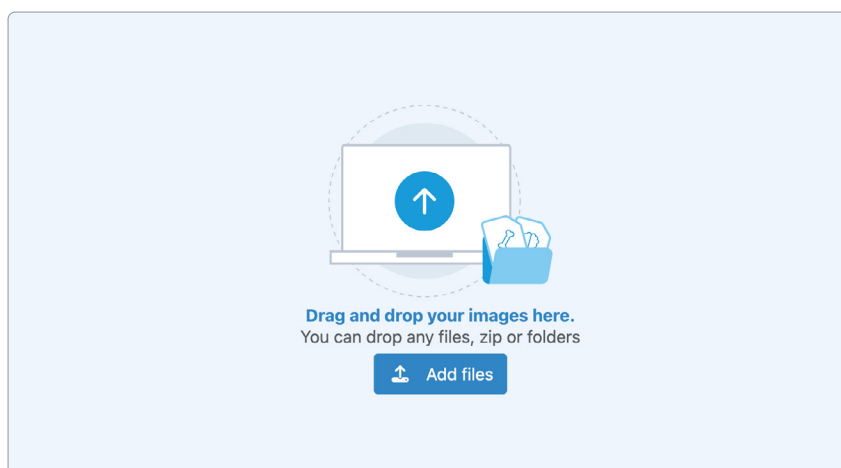

5

After completing all fields, upload DICOM images by dragging the image files to the “Upload Images” field or selecting your file location using the “Add files” button.

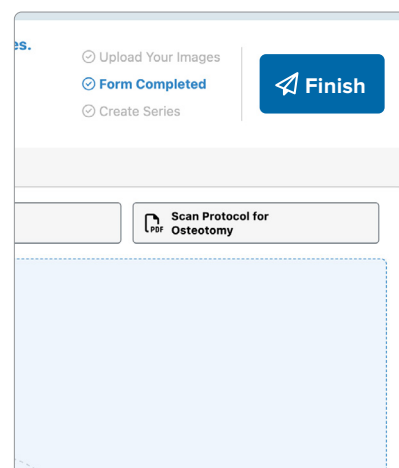

6

After uploading files, click the “Finish” button at the top right of the screen.

Note: A confirmation email with further instructions will be sent to the address provided upon login.

Document prepared by 3D-Side for clinical reference and commercial team use.
For regulatory compliance, refer to 3D-Cut FDA 510(k) K212237 cleared labeling. March 2026.



Arthrex manufacturer, authorized representative, and importer information (Arthrex eIFUs)



US patent information

© 2026-08 Arthrex, Inc. All rights reserved. LB1-000671-en-US_A

arthrex.com