

FiberSnare[®] Suture

Knee Product Highlight



The new larger loop FiberSnare suture was specifically engineered for convenience and ease of use when shuttling sutures during ligament reconstruction and repair procedures. The FiberSnare suture has an overall working length of 26 in and the construct consists of a #2 FiberWire[®] suture with a 2 in loop on one end while the opposite end is stiffened 12 in.

Key Features and Benefits

- **Versatile:** Can be used for either retrograde or antegrade drilling techniques in any orientation, allowing for delivery of either the looped or single-tail end, depending on technique
- **New color options:** Two new color options, black/white and blue/white striped, enhance visibility and allow for easier suture identification and management
- **Eliminates unloading:** Design eliminates the opportunity for the FiberStick[™] suture to unload when only a single limb is retrieved or passed, which requires knot tying or repassing the suture
- **Efficient approach:** New larger loop makes it easier to pull passing sutures through the loop

Ordering Information

Product Description	Item Number
#2 FiberSnare suture, 26 in length w/ 2 in closed loop (black/white) (a)	AR-7209SNT
#2 FiberSnare suture, 26 in length w/ 2 in closed loop (white/blue) (b)	AR-7209SNL

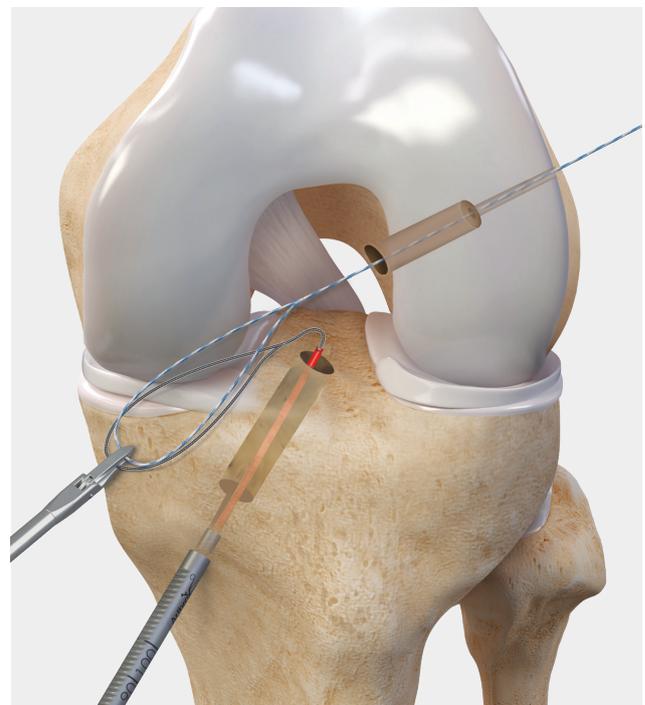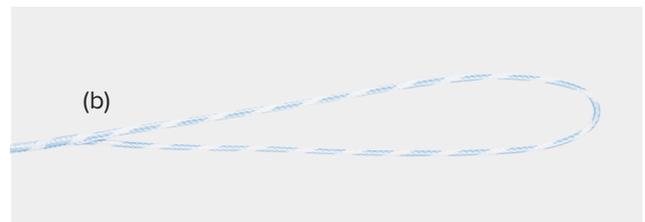

arthrex.com

© 2024-02 Arthrex, Inc. All rights reserved. DOC1-001279-en-US_A

