

DualWave™ Arthroscopy Pump

Flexibility & Versatility

Arthrex is *Raising the Bar* with the DualWave Arthroscopy Pump. The DualWave, with its direct communication with the Synergy Resection System, is an integrated inflow/outflow fluid management system that may also be used as a simple inflow-only arthroscopy pump. That means two pumps in one.

Clear, consistent turbulent-free visualization is achieved through Arthrex's latest proprietary design which precisely and independently controls the shaver or cannula outflow to maintain optimal joint distention and clarity.



Autoclavable Remote

The touch panel video display gives real time pressure and flow readings, and displays the total fluid used during a case as well as the total run time at the end of a procedure.

Flexibility is key with touch panel access and the choice of a foot pedal or autoclavable remote control that allows the user to control the DualWave from the sterile field.

User-selected defaults encompass pressure settings (by joint space), flow rate and shaver suction levels.

The Rinse and Lavage modes allow the user to control all aspects of fluid management at the touch of a button.

To help control costs, the DualWave Arthroscopy Pump can utilize Arthrex's ReDeuce™ Tubing option to dramatically lower the tubing and fluid cost per patient, while reducing room turnover times.

You will clearly see the difference with an Arthrex DualWave Arthroscopy Pump.



DualWave Arthroscopy Pump shown on optional mobile cart with Synergy Resection System. Cart may be configured for storage of fluid bags, shaver blades and/or foot pedals.



The Future of Arthroscopy is in Your Hands
<http://tinyurl.com/DUALWAVE-2013>

DualWave Arthroscopy Pump and Accessories



AR-6480 DualWave Arthroscopy Pump



AR-6482
Autoclavable Remote Control



AR-6483
Foot Pedal



AR-6481 Pump Cart

Pump Tubing




AR-6410 Main Pump Tubing



AR-6220 Extension Tubing



AR-6411 ReDeuce Pump Tubing

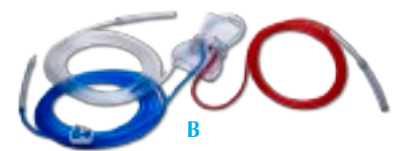

A AR-6421 ReDeuce Patient Tubing

B AR-6430 Outflow Cassette

A+B AR-6435 Outflow Tube, Set
w/ReDeuce Tubing System



AR-6215 Y-Tubing Adapter

