



ARTHREX
Graft Tubes

Clearly Better Graft Sizing and Preparation




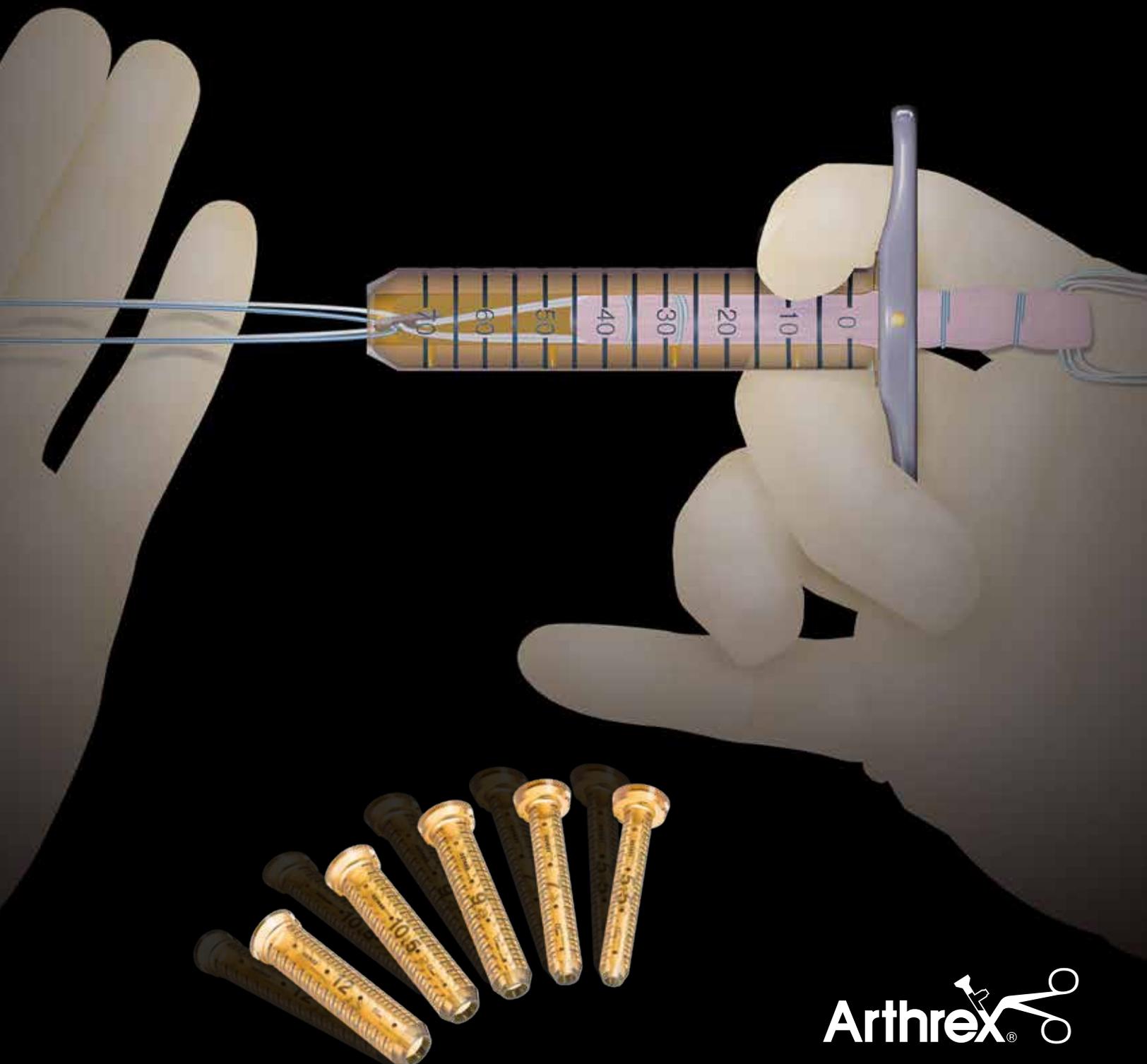


Arthrex Graft Tubes are full circumference, full length, clear tubes that facilitate graft compression, sizing and preparation. The unique transparent tubes with an etched ruler allow visualization of the graft during diameter and length sizing. Funneled entrance and attachable handle ease the entry of grafts into the sizer, allowing compression of up to 2 mm. Small holes in the Graft Tubes allow hydration of the graft or injection of biologics along the entire length. The tapered tip allows delivery of the graft directly into the tibial tunnel or medial portal.

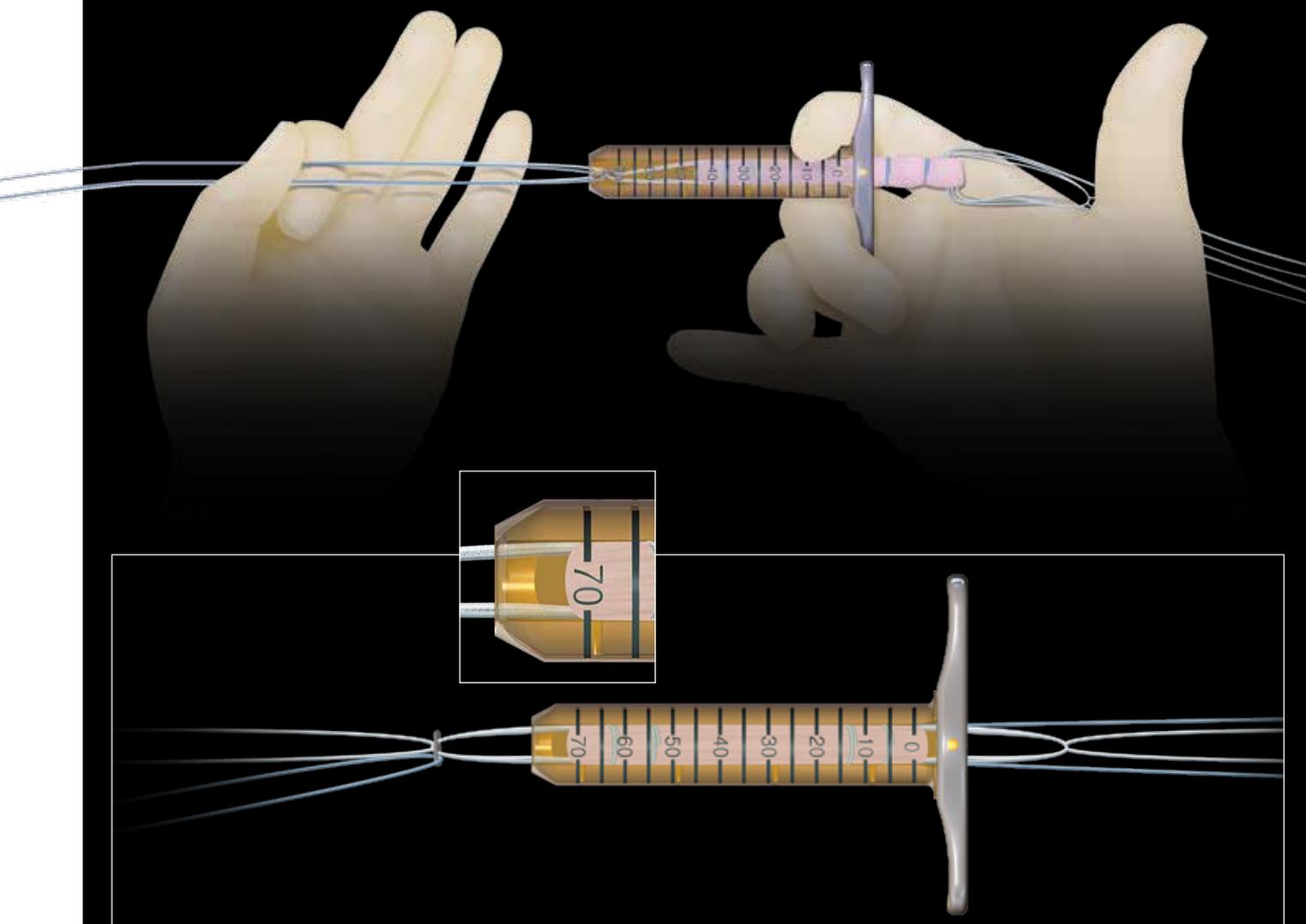


A unique transparent material with an etched ruler allow visualization of the graft during diameter and length sizing.

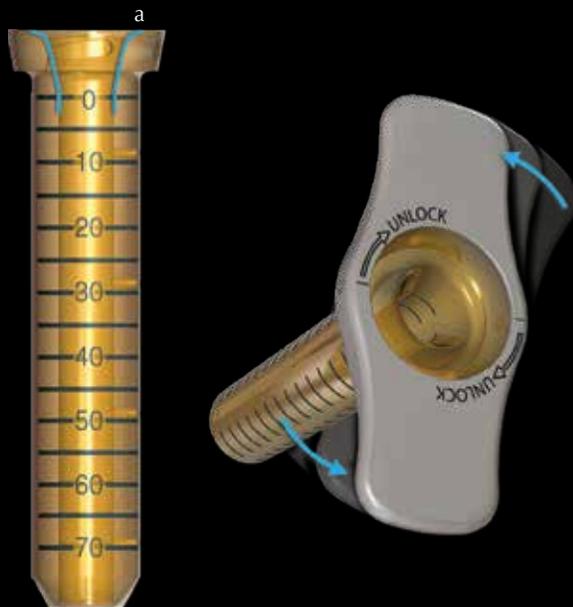

Funneled entrance (a) and attachable handle ease entry of grafts into the sizer. Grafts can be compressed up to 2 mm in diameter to facilitate passing and standardize socket sizing.

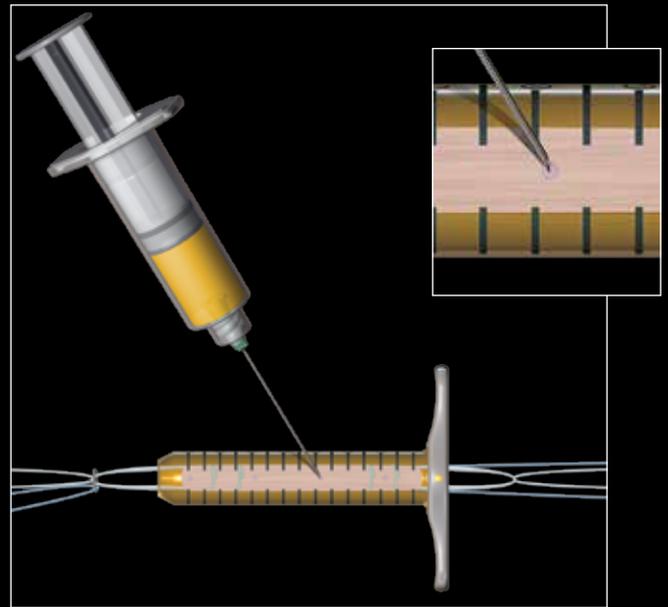

Perforations along the length of the Graft Tubes facilitate graft hydration or injection of biologics.

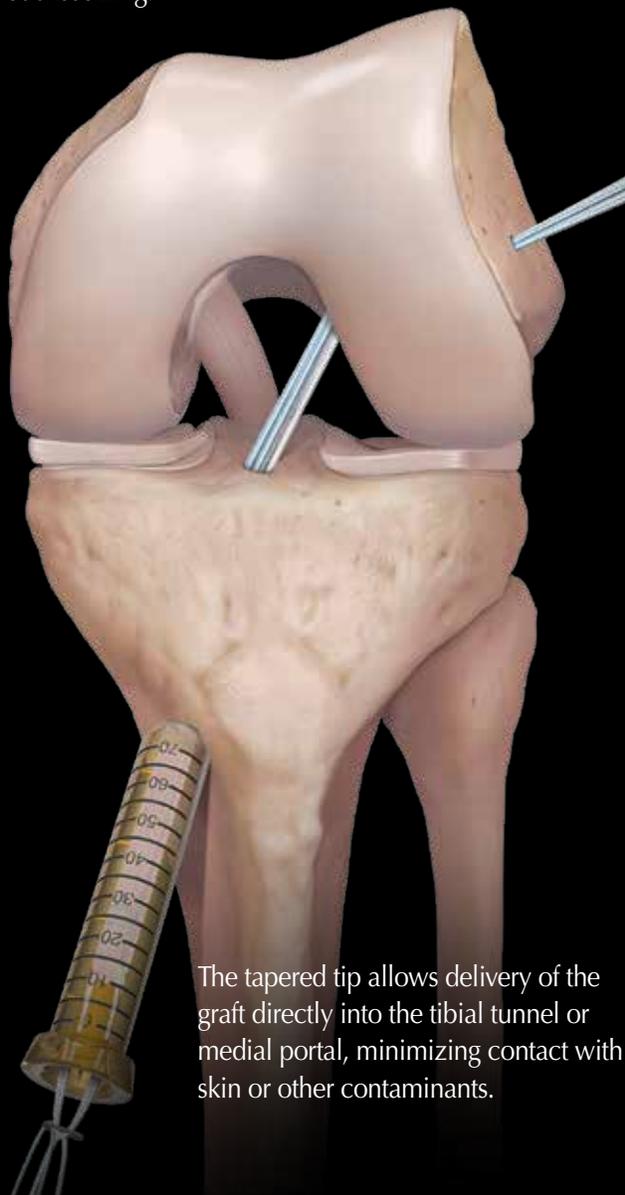

The tapered tip allows delivery of the graft directly into the tibial tunnel or medial portal, minimizing contact with skin or other contaminants.

The Graft Tube Set comes in sizes 6 mm – 13 mm in diameter including half sizes. The low profile sterilization tray can be processed independently or placed inside the RetroConstruction™ Drill Guide Instrument Set.

Ordering Information

Graft Tube Set	AR-1886-S
Graft Tubes, 5 mm – 13 mm	AR-1886-050 – 130
Graft Tube Flange	AR-1886-001
Graft Tube Instrumentation Case	AR-1886C-01



View U.S. patent information at www.arthrex.com/corporate/virtual-patent-marking

©2016, Arthrex Inc. All rights reserved. LB1-00011-EN_C