


FiberTag[®] TightRope[®] Implant

Product Highlight

The **FiberTag TightRope implant** facilitates attachment of single-ended grafts, such as quad tendon grafts, to the ACL TightRope RT and ABS implants. The FiberTag suture acts as a ripstop when placing each needle pass with the FiberLoop[®] suture, incorporating both the graft and suture.

This new implant integrates FiberTag suture with the TightRope implant for a strong, consistent connection between the suture and TightRope loop.

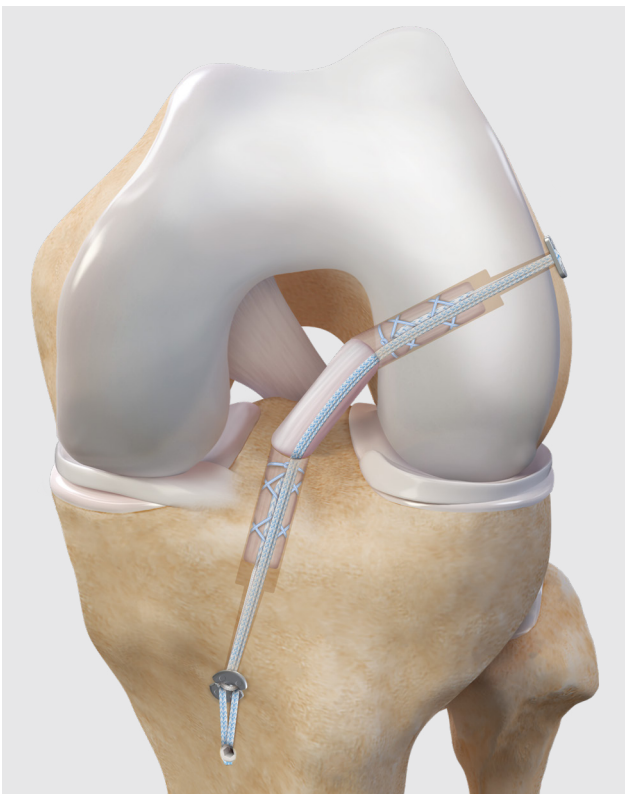

QuadLink[™] quad tendon (QT) ACL autograft with FiberTag TightRope implants and FiberTape[®] suture for *Internal/Brace*[™] technique

Features and Benefits

Improved Biomechanics


- The FiberTag rip-stop technology improves construct strength compared to conventional stitching techniques.¹ The ACL TightRope allows for precise, incremental tensioning and retensioning of the graft. Cortical fixation eliminates the risk of graft laceration from an interference screw.

Simplified Graft Preparation

- A simplified suturing technique, innovative packaging, and the new GraftClamp graft preparation instrument make preparing quadriceps tendon grafts faster and more reproducible than ever.

Optimal Healing Environment

- The new FiberTag TightRope technique reduces the overall volume of suture in the graft, creating more graft-to-bone interface and allowing the graft to circumferentially heal within the socket.

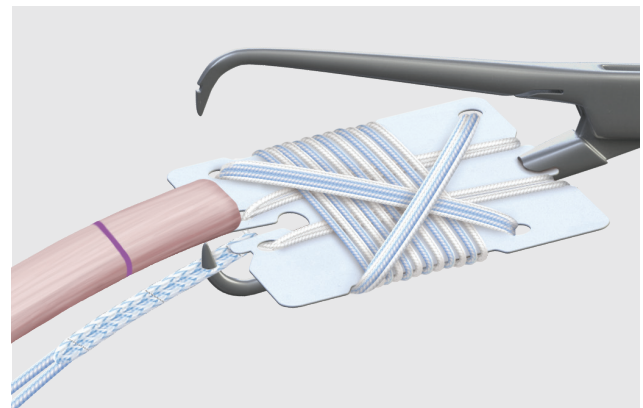

Arthrex[®]

FiberTag[®] TightRope[®] Implant

Ordering Information

FiberTag TightRope Implants

Product Description	Item Number
FiberTag TightRope and FlipCutter [®] III Implant System	AR-1288RTT-FC3
ACL FiberTag TightRope Implant	AR-1588RTT
ACL FiberTag TightRope ABS Implant	AR-1588TNT

Additional ACL Reconstruction Products

Product Description	Item Number
QuadPro [™] Tendon Harvester, 8 mm	AR-2386-08
QuadPro Tendon Harvester, 9 mm	AR-2386-09
QuadPro Tendon Harvester, 10 mm	AR-2386-10
QuadPro Tendon Harvester, 11 mm	AR-2386-11

Concave TightRope ABS Button 11 mm, round	AR-1588TB-3
Concave TightRope ABS Button 14 mm, round	AR-1588TB-4
Concave TightRope ABS Button 20 mm, round	AR-1588TB-5
ACL Backup Fixation Kit w/ BioComposite SwiveLock [®] Anchor	AR-1593-BC
ACL Backup Fixation Kit w/ PEEK SwiveLock Anchor	AR-1593-P
FlipCutter III Drill	AR-1204FF
ACL ToolBox Instrument Set	AR-1900S
PCL ToolBox Set	AR-1269S
RetroConstruction [™] Drill Guide System Instrument Set	AR-1510S
GraftClamp	AR-2386T

Products may not be available in all markets because product availability is subject to the regulatory approvals and medical practices in individual markets. Please contact your Arthrex representative if you have questions about the availability of products in your area.

Reference

1. Arthrex, Inc. Data on file (LA1-000023-en-US_A). Naples, FL; 2020.

InternalBrace[™] surgical technique is intended only to support the primary repair and is not intended as a replacement for the standard of care using biologic augmentation in a primary repair. *InternalBrace* surgical technique is intended only for soft-tissue-to-bone fixation and is not cleared for bone-to-bone fixation.